

Unit to Study a Water Meter Cutaway Model



www.edibon.com ▶PRODUCTS \$8.- FLUID MECHANICS



INTRODUCTION

The cutaway units are actual models of elements of study that have been cut to observe their interior. It allows the student to recognize all the components and understand the construction and operation of the parts that make up the element.

Water meters are usually used in residential or industrial facilities so that utility companies make the corresponding charges to users.

The fundamental difference with the gas gauges is that the meters are usually owned by the supplier and verified to ensure the correct measurement; while the gauges are usually from individuals or companies.

The water meters allow not only to know the water consumption, they also allow to evaluate if there are options for improvement or saving. It is possible to detect the existence of leaks in the network, compare our consumption with that of the average per inhabitant or measure the efficiency of the savings measures carried out as flow regulator or double-flush in the water cistern.

GENERAL DESCRIPTION

The Unit to Study a Water Meter Cutaway Model, "WMCM", allows the students to visualize the inner parts of a water meter, to evaluate the relation with each other accurately and to understand their function.

It is a model of a device that measures, registers and indicates the volume of water in cubic meters that passes through it.

It is composed of real industrial elements designed in such a way that its inner parts can be observed by the user, maintaining its mechanical and functional capacity.

1

The unit is attached to a base plate with handles for easy transport.









SPECIFICATIONS

Bench-top unit with adjustable legs.

Anodized aluminum frame and panels made of painted steel.

Fixed equipment on a base plate with handles for easy transport.

Water meter:

Material: cast iron/steel Nominal diameter: DN50. Pressure: MAP 16 bar.

Maximum working temperature: 50 °C.

Nominal flow: 250 I/min.

Manuals: This unit is supplied with the following manuals: Required Services, Assembly and Installation, Starting-up, Safety, Maintenance

& Practices Manuals.

EXERCISES AND PRACTICAL POSSIBILITIES

1.- Identification of water meters in industrial diagrams.

5.- Study and understanding of the maintenance of a water meter.

- 2.- Identification of the industrial components in an industrial plane (joint, exploded and exploded plane).
- 3.- Identification of the main components of a water meter.
- 4.- Operation of a water meter.

DIMENSIONS AND WEIGHTS

2

WMCM:

-Dimensions: 350 x 300 x 400 mm approx.


(13.77 x 11.81 x 15.74 inches approx.)

-Weight: 5 Kg approx.

(11 pounds approx.)

www.edibon.com

WMCM/ICAI. Interactive Computer Aided Instruction Software System:

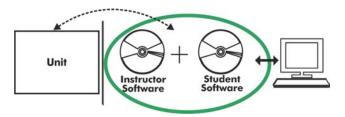

With no physical connection between unit and computer, this complete software package consists of an Instructor Software (EDIBON Classroom Manager -ECM-SOF) totally integrated with the Student Software (EDIBON Student Labsoft -ESL-SOF). Both are interconnected so that the teacher knows at any moment what is the theoretical and practical knowledge of the students.


Instructor Software

- ECM-SOF. EDIBON Classroom Manager (Instructor Software).

ECM-SOF is the application that allows the Instructor to register students, manage and assign tasks for workgroups, create own content to carry out Practical Exercises, choose one of the evaluation methods to check the Student knowledge and monitor the progression related to the planned tasks for individual students, workgroups, units, etc... so the teacher can know in real time the level of understanding of any student in the classroom.

Innovative features:

- User Data Base Management.
- Administration and assignment of Workgroup, Task and Training sessions.
- Creation and Integration of Practical Exercises and Multimedia Resources.
- Custom Design of Evaluation Methods.
- Creation and assignment of Formulas & Equations.
- Equation System Solver Engine.
- Updatable Contents.
- Report generation, User Progression Monitoring and Statistics.

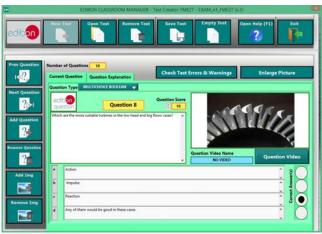

ETTE. EDIBON Training Test & Exam Program Package - Main Screen with Numeric Result Question



ECM-SOF. EDIBON Classroom Manager (Instructor Software)
Application Main Screen



ECAL, EDIBON Calculations Program Package - Formula Editor Screen



ERS. EDIBON Results & Statistics Program Package - Student Scores Histogram

Student Software

- ESL-SOF. EDIBON Student Labsoft (Student Software).

ESL-SOF is the application addressed to the Students that helps them to understand theoretical concepts by means of practical exercises and to prove their knowledge and progression by performing tests and calculations in addition to Multimedia Resources. Default planned tasks and an Open workgroup are provided by EDIBON to allow the students start working from the first session. Reports and statistics are available to know their progression at any time, as well as explanations for every exercise to reinforce the theoretically acquired technical knowledge.

Innovative features:

- Student Log-In & Self-Registration.
- · Existing Tasks checking & Monitoring.
- Default contents & scheduled tasks available to be used from the first session.
- Practical Exercises accomplishment by following the Manual provided by EDIBON.
- Evaluation Methods to prove your knowledge and progression.
- Test self-correction.
- Calculations computing and plotting.
- Equation System Solver Engine.
- User Monitoring Learning & Printable Reports.
- Multimedia-Supported auxiliary resources.


For more information see ICAI catalogue. Click on the following link: www.edibon.com/en/files/expansion/ICAI/catalog




ERS. EDIBON Results & Statistics Program Package - Question Explanation



ESL-SOF. EDIBON Student LabSoft (Student Software)
Application Main Screen

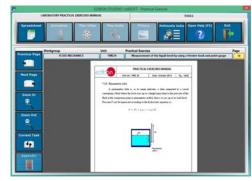

EPE. EDIBON Practical Exercise Program Package Main Screen

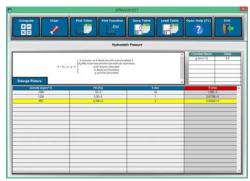

ECAL. EDIBON Calculations Program Package Main Screen

* Specifications subject to change without previous notice, due to the convenience of improvement of the product.



C/ Julio Cervera, 10-12-14. Móstoles Tecnológico. 28935 MÓSTOLES. (Madrid). ESPAÑA - SPAIN. Tel.: 34-91-6199363 Fax: 34-91-6198647

E-mail: edibon@edibon.com Web: www.edibon.com

Edition: ED01/18 Date: August/2018

